

PREGNANCY TEST

What makes

True[®] 20 & True[®] 20 Plus
Truly Exceptional?



Exceptional Quality
Exceptional Price

CLIA '88
WAIVED

True[®] 20 & True[®] 20 Plus

*Truly Accurate, Truly Economical,
Truly Exceptional*



The World's Source for Clinical Diagnostics

PREGNANCY TEST

EVERYTHING YOU NEED AT A PRICE YOU CAN'T RESIST

Everything you need is included in one foil pouch (including dropper). The device features a clear, easy-to-read background that ensures reliable results each and every time, providing the highest quality at an affordable price. You will *truly* notice the difference.



True[®] 20 & True[®] 20 Plus Benefits

Enhanced Readability - Background clears within seconds

Sensitive - 20 mIU/mL hCG for BOTH urine and serum samples

Fast - Positive results can often be seen in as little as 10 seconds

Small Sample Size - Only 2 drops of sample (urine or serum) required

Reliable, Specific, and Accurate - Specific mono and polyclonal antibodies ensure early pregnancy detection while preventing cross reactivity with other hormones present in patient samples

Convenient - Room temperature storage

CLIA Waived status for urine samples (CPT Code: 81025QW)

Products

True[®] 20 hCG (Urine)

True[®] 20 Plus hCG (Serum/Urine)

Size

50 Tests

50 Tests

Reference No.

1430-050

1440-050

General Test Procedure


1. Remove device and disposable dropper from foil pouch.
2. Label the test device with patient or control identification.
3. Dispense 2 drops of the specimen into the round sample well.
4. Read result as positive or negative.



Positive



Negative

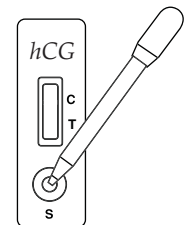

The World's Source for Clinical Diagnostics

- General & Special Chemistry
- Pregnancy
- Point-of-Care
- Rapid Tests
- Therapeutic Drug Monitoring
- Serology
- Occult Blood
- Urinalysis
- Clinical Instrumentation
- Pipettors

1261 North Main Street
Boerne, Texas 78006 USA

Toll Free: 866.STANBIO
(782.6246)

Tel: 830.249.0772
Fax: 830.249.0851

stanbio@stanbio.com
www.stanbio.com



DISTRIBUTED BY

Lit # RSS0012.00
©2007 Stanbio Laboratory
Printed in USA.