

PREGNANCY TEST



Save Green:

First-ever Earth Friendly & Budget Friendly pregnancy test

CLIA '88
WAIVED

alpha[®]

*Less Packaging, Less Waste,
More Value*



STANBIO
LABORATORY

The World's Source for Clinical Diagnostics

PREGNANCY TEST

EVERYTHING YOU NEED AND MORE THAN YOU HOPED FOR

When we say eco-friendly, we mean it as a double threat. Not just economically friendly, but also ecologically friendly - which is a hard balance to strike.

First, alpha® is gentle on your budget. By bundling our pouched tests in quantities of 100 in a resealable bag, we save on costly custom-printed packaging. Individually wrapped dipstrips prolong the life of the tests, saving you wasted product and money.

We also aspired to be gentle to the environment. The resealable bags that alpha® comes packaged in are fully recyclable and/or reusable around your office or lab. The cassettes undergo less processing, which saves energy and by-products. Even this very sheet of paper was printed on 100% biodegradable stone-based paper, which doesn't require any chemical treatments, bleaching, or washing - all of which saves trees, water, chemical runoff, and energy. It even uses less ink during the printing process.

Stanbio Laboratory has always been proud to offer the very best diagnostic products on the market in order to improve the health of patients. Now we're raising the bar and minimizing our environmental impact to improve the health of our planet, too.



alpha® Benefits

Enhanced Readability - Background clears within seconds

Sensitive - 25 mIU/mL hCG for BOTH urine and serum samples

Fast - Positive results can often be seen in as little as 10 seconds

Small Sample Size - Only 4 drops of sample (urine or serum) required



Reliable, Specific, and Accurate - Specific mono and polyclonal antibodies ensure early pregnancy detection while preventing cross reactivity with other hormones present in patient samples

Convenient - Room temperature storage (15°C - 30°C)

CLIA Waived status for urine samples (CPT Code: 81025QW)

General Test Procedure

1. Remove device and disposable dropper from foil pouch.
2. Label the test device with patient or control identification.
3. Dispense 4 drops of the specimen into the round sample well.
4. Read result as positive or negative.

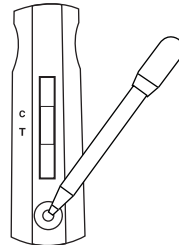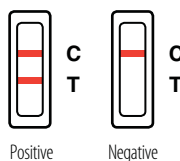

The World's Source for Clinical Diagnostics

- General & Special Chemistry
- Pregnancy
- Point-of-Care
- Rapid Tests
- Therapeutic Drug Monitoring
- Serology
- Occult Blood
- Urinalysis
- Clinical Instrumentation
- Pipettors

1261 North Main Street
Boerne, Texas 78006 USA
Toll Free: 866.STANBIO
(782.6246)
Tel: 830.249.0772
Fax: 830.249.0851
stanbio@stanbio.com
www.stanbio.com



DISTRIBUTED BY

Item # RSS0032
Lit # RSS0032.00
Eff. Date: 7/16/10
DAR: 10.1979
©2010 Stanbio Laboratory